


## Waterbath WNB 22

The combination of corrosion-resistant stainless steel, precise electronics and multiple temperature protection guarantees the highest level of safety in the laboratory.



On this page, you can find all the essential technical data on the Memmert water bath. Our customer relations team will be pleased to help if you want further information. If you should require a customised special solution, please contact our technical specialists at [sales@memmert.com](mailto:sales@memmert.com).

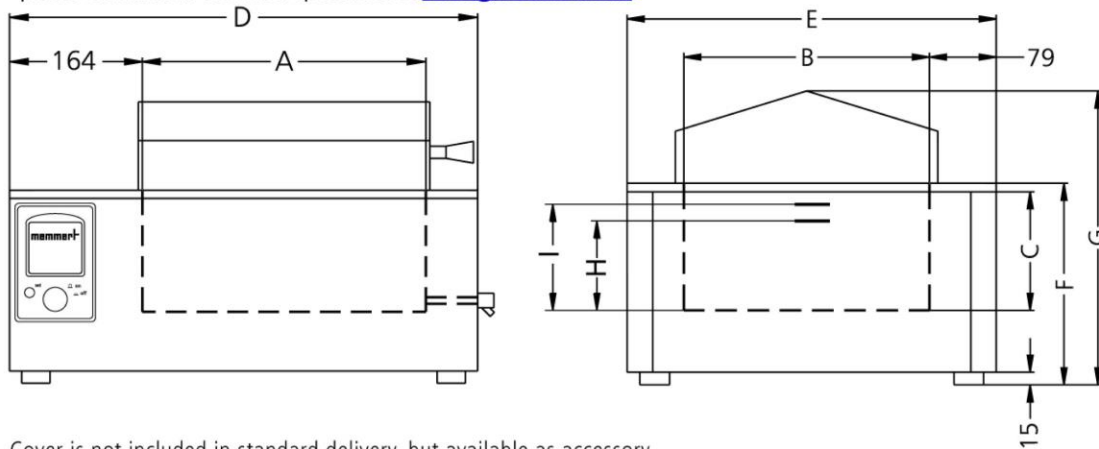

Cover is not included in standard delivery, but available as accessory.



تهران، خیابان مطهری، خیابان فجر، پلاک 38



021-88832100

021-88321727

021-88321728



[www.abzarteb.com](http://www.abzarteb.com)

[abzar.teb](https://www.instagram.com/abzar.teb)

## Temperature

Setting temperature range	+10 to +95 °C and boiling stage
Working temperature range in °C	min. 5 above ambient up to +95 °C with additional boiling mode

## Control of standard components

Temperature sensor	1 Pt100 sensor class A in 4-wire-circuit
Timer	integrated digital timer from 1 min. to 99,99 hours for:ON continuous operationWAIT (delayed on for continuous andlimited timed operation)HOLD

## Control technology

Controller	digital display of set and actual temperature and of (remaining) programme time
------------	---

## Safety

Temperature control	mechanical temperature limiter TB protectionclass 1 switching the heating off at approx. 30°Cabove max. temperature of the bath
Temperature control	in case of overtemperature due to failure,the heating is switched off at approx. 10°Cabove the set temperature (fixed value)
Autodiagnostic system	microprocessor PID-temperature controllerwith integrated autodiagnostic system withfault indicator

## Heating concept

Heating Baths	corrosion-proof large-area heating onthree sides
---------------	--

## Stainless steel interior

Dimensions W x H x D in mm	$w_{(A)} \times h_{(C)} \times d_{(B)}$ : 350 x 220 x 290 mm
Interior	easy-to-clean interior, made of stainlesssteel, reinforced by deep drawn ribbing,material no. 1.4301 (ASTM 304), laser-welded
Volume	22 l
Liquid level min.	177 mm
Liquid level max.	200 mm

## Textured stainless steel casing

Dimensions	$w_{(D)} \times h_{(G)} \times d_{(E)}$ : 578 x 405 x 436 mm
------------	--



تهران، خیابان مطهری، خیابان فجر، پلاک 38



021-88832100

021-88321727

021-88321728



www.abzarteb.com



abzar.teb

**Electrical data**

Voltage	230 V, 50/60 Hz
Electrical load	approx. 2000 W (during heating)

**Ambient conditions**

Installation	The vent openings in the left and back side must remain unobstructed. Minimum wall spacing on all sides is 80 mm. The minimum spacing from the top of the bath to the next ceiling is 750 mm.
Ambient temperature	+5 °C to +40 °C
Humidity rh	max. 80 %, non-condensing
Overvoltage category	II
Pollution degree	2

**Packing/shipping data**

Transport information	The appliances must be transported upright
Customs tariff number	8419 8998
Country of origin	Federal Republic of Germany
WEEE-Reg.-No.	DE 66812464
Dimensions approx incl. carton	w x h x d: 670 x 530 x 460 mm
Net weight	approx. 16 kg
Gross weight carton	approx. 23 kg

Standard units are safety-approved and bear the test marks



تهران، خیابان مطهری، خیابان فجر، پلاک 38



021-88832100

021-88321727

021-88321728



www.abzarteb.com



abzar.teb